

Mimetix[®] Aligned 384-well plate with rhodamine

AMS.TECL-015-1X AMS.TECL-015-8X



The Mimetix Aligned 384-well plate is an easy to use tool for the culture of cells which are influenced by topographical features, including cells that myelinate. The fibres contain rhodamine for ease of imaging and analysis.

Product Description:

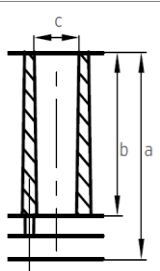
Mimetix aligned fibres containing rhodamine-6G are incorporated into a standard 384-well plate frame using a proprietary laser-welding technology which provides minimal base distortion and avoids the use of glues.

The base is 190 µm polystyrene with excellent optical properties and high light transmission.

Scaffold Specifications:

- Material: medical-grade poly-L-lactide (PLLA)
- Label: Rhodamine-6G
- Orientation: Aligned
- Fibre diameter: 2 µm
- Scaffold thickness: 2-4 µm
- Scaffold density: 130 fibres/mm
- Non-biodegradable in *in vitro* applications

Plate Specifications

	Plate width:	127.76 ± 0.25 mm
	Plate depth:	85.48 ± 0.25 mm
	Plate height (a):	14.35 ± 0.25 mm
	Well depth (b):	11.35 ± 0.25 mm
	Well diameter (c):	3.70 ± 0.10 mm
	Distance to centre of A1 from top edge:	8.99 ± 0.25 mm
	Distance to centre of A1 from left edge:	12.13 ± 0.25 mm
	Pitch (distance between A1 and A2):	4.50 mm

- Supplied with lid and individually-sealed
- Treated with gamma or e-beam irradiation
- Store at room temperature in the dark
- Manufactured in the U.K.

Features:

- Compatible with industry-standard automated handling and imaging equipment including fluorescence microscopy
- Protocols for cell seeding, assays, and imaging are available

AMSBIO | www.amsbio.com | info@amsbio.com

 **UK & Rest of the World**
AMS Biotechnology (Europe) Ltd
184 Park Drive, Milton Park
Abingdon, UK
T: +44 (0)1235 828 200
F: +44 (0) 1235 820 482

 **North America**
amsbio LLC
1035 Cambridge Street,
Cambridge, MA 02141
T: +1 (617) 945-5033 or
T: +1 (800) 987-0985
F: +1 (617) 945-8218

 **Germany**
AMS Biotechnology (Europe) Ltd
Bockenheimer Landstr. 17/19
60325 Frankfurt/Main
T: +49 (0) 69 779099
F: +49 (0) 69 13376880

 **Switzerland**
AMS Biotechnology (Europe) Ltd
Centro Nord-Sud 2E
CH-6934 Bioggio-Lugano
T: +41(0) 91 604 55 22
F: +41(0) 91 605 17 85



Mimetix[®] Aligned 384-well plate with rhodamine

Precondition

The Mimetix scaffold needs to be wetted with ethanol in order to allow the cells to attach to the fibres.

- Add 30 μ L of 20% ethanol per well.
- Allow ethanol to soak into the membrane for 5 min, then aspirate it carefully without touching the scaffold.

Wash

- Wash scaffold twice with PBS.
- Leave scaffold in cell culture medium until cell seeding.

Seed

Seeding densities are guidelines only. Density as for 2D culture.

- Add your cells suspended in 30-60 μ L cell culture medium.

Exchange medium

- For long-term experiments semi-exchange the cell culture medium every 3 days.

AMSBIO | www.amsbio.com | info@amsbio.com



UK & Rest of the World
AMS Biotechnology (Europe) Ltd
184 Park Drive, Milton Park
Abingdon, UK
T: +44 (0)1235 828 200
F: +44 (0) 1235 820 482



North America
amsbio LLC
1035 Cambridge Street,
Cambridge, MA 02141
T: +1 (617) 945-5033 or
T: +1 (800) 987-0985
F: +1 (617) 945-8218



Germany
AMS Biotechnology (Europe) Ltd
Bockenheimer Landstr. 17/19
60325 Frankfurt/Main
T: +49 (0) 69 779099
F: +49 (0) 69 13376880



Switzerland
AMS Biotechnology (Europe) Ltd
Centro Nord-Sud 2E
CH-6934 Bioggio-Lugano
T: +41(0) 91 604 55 22
F: +41(0) 91 605 17 85